# microlife®

## BP B3 Comfort PC Blood Pressure Monitor

EN → 1	CZ → 38	HU → 70
RU → 9	SK → 46	HR → 79
BG → 19	SL → 54	PL → 87
RO → 29	SR → 62	

IB BP B3 Comfort PC E-V11 4823  
Revision Date: 2023-11-06

Microlife Corporation  
9F, 431, RuiGuang Road, NeiHu  
Taipei 11492, Taiwan, China  
www.microlife.com

Microlife UAB  
P. Lukšio g. 32,  
08222 Vilnius  
Lithuania

CE0044

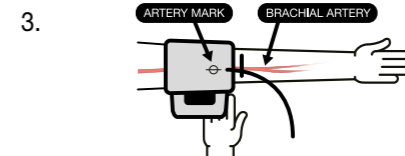
### Preparation

Sit on a back-supported chair and keep your legs uncrossed. / Сядьте, не скрещивая ноги, на стул со спинкой. / Сядьте на стол с облегалка и не кръстосвайте краката си. / Așezajți-vă într-un scaun cu spătar și nu încrucișați picioarele / Sediñete si na židli s opěrkou a nepřekřížujte nohy. / Sadnite si na stoličku s opěrkou a nepřekřížujte si nohy. / Usedite se na stol s podprtím hrbtom in nogami na tleh (ne prekrížajte jih). / Sedite na stolicu koja podupire leđa i nemojte prekrštati noge. / Uljön egy háttámlás székre és lábait ne rakja keresztbe! / Sjednite na stolicu koja podupire leđa i neka Vam noge ne budu prekrížene. / Usiądź na krzesle wspieranym plecami i nie rozstawiaj nóg.

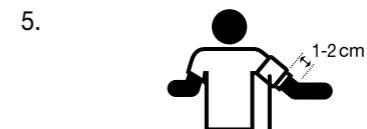
Avoid thick or close-fitting garments on the upper arm. / Освободите плечо от плотной или плотно облегающей одежды. / Избягвайте дебели или плътно прилепващи дрехи върху горната част на ръката. / Evitați îmbracaminte groasă sau stâmta pe braț / Vyhňte se těsnému nebo přiléhavému oděvu na horní části paže. / Vyhňte sa těsnému alebo priléhavému odevu na hornej časti ramena. / Izogibajte se tesnim oblačilom, ki stisnejo nadlaket. / Izbegavajte usku odeću na nadlaktici. / Kerülje a vastag vagy szoros ruhákat a felkaron! / Izbjegavajte usku odeću na nadlaktici. / Unikaj grubej lub dopasowanej odzieży na ramieniu.

Place the artery-mark on the cuff over your artery. / Поместите манжету так, чтобы значок артерии на манжете находился над артерией руки в сгибе локтя. / Поставете маркера на маншета върху артерията. / Plasați semnul artera al manșetei peste artera dumneavoastră / Umístěte manžetu se značkou tepny nad Vaši tepnu. / Namestite znak za arterijo, ki je na manšeti, na vašo arterijo. / Postavite oznaku za arteriju na manžetni preko Vaše arterije. / Helyezze a mandzsetta artériajelzójét az artériája fölé! / Postavite oznaku arterije na manžeti preko Vaše arterije. / Umieść mankiecie na tętnicy.

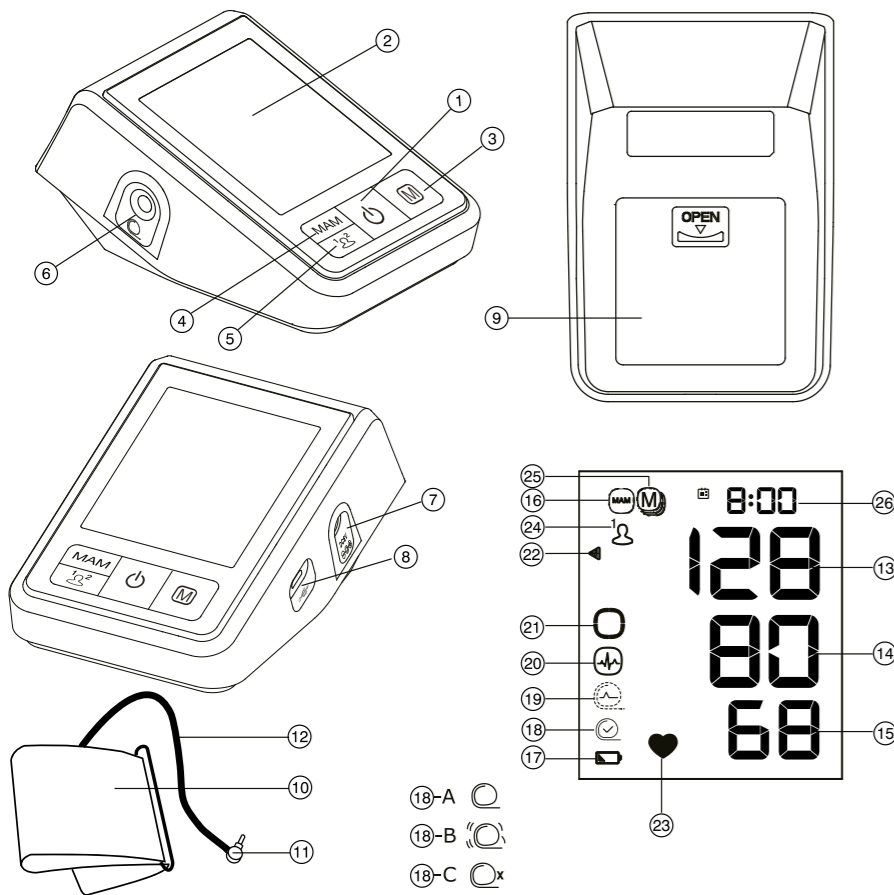
Fit the cuff closely, but not too tight. / Поместите манжету так, чтобы она плотно, но не туго прилежала. / Стегнете добре маншета, но не твърде много. / Fixați manșeta ferm dar nu prea strânsă / Nasade manžetu těsně, ale ne příliš těsně. / Nasade manžetu tesne, ale nie příliš tesne. / Manšeto dobro namestite, vendar ne pretisno. / Dobro zategnite manžetnu, ali ne previše stegnuto. / Rögzítse a mandzsettát szorosan, de ne feszítse meg túl erősen! / Postavite manžetu blizu, ali ne pretisno. / Zaciñnij mankiet dokladnie, jednak niezbyt silnie.

Position the cuff 1-2 cm above your elbow. / Расположите манжету на 1-2 см выше локтя. / Поставете маншета 1-2 см над лакътя. / Plasați manșeta cu 1-2 cm deasupra articulației / Manžetu umístěte 1-2 cm nad loket. / Manžetu umestnite 1-2 cm nad laket. / Namestite manšeto 1-2 cm nad komolcem. / Postavite manžetnu 1-2 cm iznad Vašeg lakta. / A mandzsetta 1-2 centiméterrel legyen a könyökétől feljebb! / Postavite manžetu 1-2 cm iznad Vašeg lakta. / Ustaw mankiet na 1-2 cm powyżej łokcia.

Keep your arm still and do not speak during the measurement. / Держите руку неподвижно и не разговаривайте во время измерения. / Дръжте ръката си неподвижна и не говорете по време на измерването. / Nu mișcați mâna și nu vorbiți în timpul măsurării tensiunii arteriale. / Držte ruku uvolněnou a během měření nemlve. / Držte ruku uvolněnou a počas merania nerozprávjajte. / Med merjenjem ne premikajte roke in ne govorite. / Držite ruku mimo i nemojte pričati tokom merenja. / Ne mozgassa a karját és ne beszéljen mérés közben! / Držite ruku mimo i nemojte pričati tijekom mjerenja. / Trzymaj rękę nieruchomo i nie mów podczas pomiaru.

1.

Avoid eating, bathing, smoking or caffeine (approx. 30 min).  
 Избегайте приема пищи, купания, курения или приема кофеиносодержащих продуктов (прибл. за 30 минут до измерения).  
 Избягвайте ядене, къпане, пушене или кофеин (около 30 минути).  
 Evitați sa mâncați, sa faceți baie, sa fumați, sau sa beți cafea (aproximativ 30 min).  
 Vyhněte se jídlu, koupání, kouření nebo kofeinu (přibližně 30 minut).  
 Vyhňte sa jedlu, kúpaniu, fajčeniu alebo kofeínu (približne 30 minút).  
 Izogibajte se hrani, kupanju in kajeju (približno 30 minut).  
 Izbegavajte hranu, kupanje, pušenje i kofein (približno 30 minuta).  
 Kerülje az evést, fürdést, dohányzást vagy koffeinfogyasztást (a megelőző kb. 30 percben)!  
 Izbjegavajte jesti, kupati se, pušiti i kofein (oko 30 min.)  
 Unikaj jedzenia, kąpieli, palenia lub kofeiny (około 30 minut).

2.

Avoid activity and relax for 5-10 min.  
 В течение 5-10 мин. оставайтесь в расслабленном состоянии и избегайте физической нагрузки.  
 Избягвайте активности и се отпуснете в продължение на 5-10 минути.  
 Evitați orice activitate și relaxați-vă 5-10 min.  
 Vyhněte se fyzické aktivitě a odpočívajte 5-10 minut.  
 Vyhňte sa fyzickej aktivite a odpočívajte 5-10 minút.  
 Izogibajte se fizični aktivnosti in počivajte 5-10 minut.  
 Izbegavajte fizičku aktivnost i opustite se tokom 5-10 minuta.  
 Kerülje a fizikai tevékenységet, pihenjen 5-10 percig!  
 Izbjegavajte aktivnost i opustite se 5-10 min.  
 Unikaj aktywności i zrelaksuj się przez 5-10 minut.

3.

Measure before medication intake.  
 Измеряйте до приема лекарства.  
 Измерете преди приема на лекарството ви.  
 Măsurati înainte de a lua medicația.  
 Měření proveďte před přijetím léku.  
 Meranie uskutočnite pred užitím liekov.  
 Merjenje izvedite pred zaužitjem zdravil.  
 Obavite merenje pre uzimanja leka.  
 Mérjen a gyógyszer bevétele előtt!  
 Izvedite mjerjenje prije nego što uzmete lijek.  
 Zmierz przed przyjęciem leku.

Name of Purchaser / Ф.И.О. покупателя /  
 Име на купувача / Numele cumpărătorului /  
 Jméno kupujícího / Meno zákazníka / Ime in  
 priimek kupca / Ime i prezime kupca / Vásárló  
 neve / Ime i prezime kupca / Imię i nazwisko  
 nabywcy

Serial Number / Серийный номер / Серийен  
 номер / Număr de serie / Výrobní číslo /  
 Výrobné číslo / Serijska številka / Serijski broj /  
 Sorozatszám / Serijski broj / Numer seryjny

Date of Purchase / Дата покупки / Дата на  
 закупуване / Data cumpărării / Datum nákupu /  
 Dátum kúpy / Datum nakupa / Datum kupovine /  
 Vásárlás dátuma / Datum kupovine / Data zakupu

Specialist Dealer / Специализированный дилер /  
 Специалист дистрибутор / Distribuitor de spe-  
 cialitate / Specializovaný dealer / Špecializovaný  
 predajca / Specjalizirani trgovac / Ovlašćeni diler /  
 Forgalmazó / Ovlašteni prodavač / Przedstawiciel

- ① ON/OFF button
- ② Display
- ③ M-button (memory)
- ④ MAM button
- ⑤ User Button
- ⑥ Cuff socket
- ⑦ Mains Adapter Socket
- ⑧ USB Port
- ⑨ Battery compartment
- ⑩ Cuff
- ⑪ Cuff connector
- ⑫ Cuff tube

## Display

- ⑬ Systolic value
- ⑭ Diastolic value
- ⑮ Pulse rate
- ⑯ MAM Mode
- ⑰ Battery display
- ⑱ Cuff fit check
  - A: Suboptimal cuff fit
  - B: Arm movement indicator «**Err 2**»
  - C: Cuff pressure check «**Err 3**»
- ⑲ Cuff signal indicator «**Err 1**»
- ⑳ Irregular heartbeat (IHB) symbol
- ㉑ MAM Interval time
- ㉒ Traffic light indicator
- ㉓ Pulse indicator
- ㉔ User indicator
- ㉕ Stored value
- ㉖ Date/Time

Read the important information in these instructions for use before using this device. Follow the instructions for use for your safety and keep it for future reference.

Type BF applied part

Keep dry

Manufacturer

Batteries and electronic devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

Authorized representative in the European Community

Catalogue number

Serial number  
(YYYY-MM-DD-SSSSS;  
year-month-day-serial number)

Caution

Humidity limitation for operating **and** storage

Temperature limitation for operating **or** storage

Medical device

Keep away from children of age 0 - 3

CE Marking of Conformity

## Intended use:

This oscillometric blood pressure monitor is intended for measuring non-invasive blood pressure in people aged 12 years or older.

Dear Customer,

This device was developed in collaboration with physicians and clinical tests carried out prove its measurement accuracy to be of a very high standard.\*

If you have any questions, problems or want to order spare parts please contact your local Microlife-Customer Service. Your dealer or pharmacy will be able to give you the address of the Microlife dealer in your country. Alternatively, visit the internet at [www.microlife.com](http://www.microlife.com) where you will find a wealth of invaluable information on our products.

Stay healthy – Microlife Corporation!

\* This device is tested according to the ESH protocol and ISO81060-2:2013.

## Table of contents

### 1. Using the device for the first time

- Inserting the batteries
- Setting the date and time
- Selecting the correct cuff
- Selecting the user
- Selecting standard or MAM mode

### 2. Checklist for taking a reliable measurement

### 3. Taking a blood pressure measurement

- How not to store a reading
- How do I evaluate my blood pressure
- Appearance of the irregular heartbeat (IHB) symbol

### 4. Data memory

- Viewing the stored values
- Clearing all values

### 5. Battery indicator and battery change

- Low battery
- Flat battery - replacement
- Which batteries and which procedure?
- Using rechargeable batteries

### 6. Using a mains adapter

### 7. PC-Link functions

### 8. Error Messages

### 9. Safety, care, accuracy test and disposal

- Device care
- Cleaning the cuff
- Accuracy test
- Disposal

### 10. Guarantee

### 11. Technical Specifications

## 1. Using the device for the first time

### Inserting the batteries

After you have unpacked your device, first insert the batteries. The battery compartment ⑨ is on the bottom of the device. Insert the batteries (4 x 1.5 V, size AA), thereby observing the indicated polarity.

### Setting the date and time

- After the new batteries are fitted, the year number flashes in the display. You can set the year by pressing the M-button ③. To confirm and then set the month, press the MAM button ④.
- Press the M-button to set the month. Press the MAM button to confirm and then set the day.
- Follow the instructions above to set the day, hour and minutes.
- Once you have set the minutes and pressed the MAM button, the date and time are set and the time is displayed.
- If you want to change the date and time, press and hold the MAM button for approx. 3 seconds until the year number starts to flash. Now you can enter the new values as described above.

### Selecting the correct cuff

Microlife offers different cuff sizes. Select the cuff size to match the circumference of your upper arms (measured by close fitting in the centre of the upper arm).

Cuff size	for circumference of upper arm
S	17 - 22 cm
M	22 - 32 cm
M - L	22 - 42 cm
L	32 - 42 cm
L - XL	32 - 52 cm

 Pre-shaped cuffs are optionally available.

 Only use Microlife cuffs.

- ▶ Contact your local Microlife Service if the enclosed cuff ⑩ does not fit.
- ▶ Connect the cuff to the device by inserting the cuff connector ⑪ into the cuff socket ⑥ as far as it will go.

 If you buy a spare Microlife cuff, please remove the cuff connector ⑪ from the cuff tube ⑫ from the cuff supplied with the original device and insert this cuff connector into the tube of the spare cuff (valid for all cuff sizes).

## Selecting the user

This device allows to store the results for 2 individual users.

- ▶ Select the intended user (user 1 or user 2 (24)) by pressing the user button (5).

 Before each measurement, ensure that the correct user is selected.

## Selecting standard or MAM mode

Before each measurement, select standard (single measurement) or MAM mode (automatic triple measurement). In MAM mode, 3 measurements are automatically taken in succession and the result is then automatically analysed and displayed. Because the blood pressure constantly fluctuates, a result obtained in this way is more reliable than when a single measurement is performed.

- To select MAM mode, press the MAM button (4) until the MAM-symbol (16) appears on the display. To change to standard mode (single measurement), press the MAM-button again, until the MAM-symbol disappears.
- The bottom, right hand section of the display shows a 1, 2 or 3 to indicate which of the 3 measurements is currently being taken.
- There is a break of 15 seconds between the measurements. A count down indicates the remaining time.
- The individual results are not displayed. Your blood pressure will only be displayed after all 3 measurements are taken.
- Do not remove the cuff between measurements.
- If one of the individual measurements was questionable, a fourth one is automatically taken.

## 2. Checklist for taking a reliable measurement

- ▶ Avoid activity, eating or smoking immediately before the measurement.
- ▶ Sit down on a back-supported chair and relax for 5 minutes. Keep your feet flat on the floor and do not cross your legs.
- ▶ **Always measure on the same arm** (normally left). It is recommended that doctors perform double arm measurements on a patients first visit in order to determine which arm to measure in the future. The arm with the higher blood pressure should be measured.
- ▶ Remove close-fitting garments from the upper arm. To avoid constriction, shirt sleeves should not be rolled up - they do not interfere with the cuff if they are laid flat.
- ▶ Always ensure that the correct cuff size is used (marking on the cuff).

- Fit the cuff closely, but not too tight.
- Make sure that the cuff is positioned 1-2 cm above the elbow.
- The **artery mark** on the cuff (ca.3 cm long bar) must lie over the artery which runs down the inner side of the arm.
- Support your arm so it is relaxed.
- Ensure that the cuff is at the same height as your heart.

## 3. Taking a blood pressure measurement

1. Select standard (single measurement) or MAM mode (automatic triple measurement): see details in chapter 1.»
2. Press the ON/OFF button (1) to start the measurement.
3. The cuff will now pump up automatically. Relax, do not move and do not tense your arm muscles until the measurement result is displayed. Breathe normally and do not talk.
4. The cuff fit check (18) on the display indicates that the cuff is perfectly placed. If the icon (18)-A appears, the cuff is fitted suboptimally, but it is still ok to measure.
5. The measurement is performed during the inflation. The inflation speed may vary, this is a normal occurrence.
6. During the measurement, the pulse indicator (23) flashes in the display.
7. The result, comprising the systolic (13) and the diastolic (14) blood pressure and the pulse rate (15) are displayed. Note also the explanations on further display symbols in this booklet.
8. When the device has finished measuring, remove the cuff.
9. Switch off the device. (The monitor does switch off automatically after approx. 1 min.)

 You can stop the measurement at any time by pressing the ON/OFF button (e.g. if you feel uneasy or an unpleasant pressure sensation).

### How not to store a reading

As soon as the reading is displayed press and hold the ON/OFF button (1) until «M» (25) is flashing. Confirm to delete the reading by pressing the MAM button (4).

 «CL» is displayed when the reading is deleted from the memory successfully.

### How do I evaluate my blood pressure

The triangle on the left-hand edge of the display (22) points at the range within which the measured blood pressure value lies. The value is either within the optimum (green), elevated (yellow) or high (red) range.

The classification corresponds to the following ranges defined by international guidelines (ESH, ESC, JSH). Data in mmHg.

Range	Systolic	Diastolic	Recommendation
1. blood pressure too high	≥135	≥85	Seek medical advice
2. blood pressure elevated	130 - 134	80 - 84	Self-check
3. blood pressure normal	<130	<80	Self-check

The higher value is the one that determines the evaluation.

Example: a blood pressure value of **140/80** mmHg or a value of **130/90** mmHg indicates «blood pressure too high».

### Appearance of the irregular heartbeat (IHB) symbol

This symbol  indicates that an irregular heartbeat was detected. In this case, the measured blood pressure may deviate from your actual blood pressure values. It is recommended to repeat the measurement.

### Information for the doctor in case of repeated appearance of the IHB symbol:

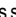
This device is an oscillometric blood pressure monitor that also measures the pulse during blood pressure measurement and indicates when the heart rate is irregular.

## 4. Data memory

This device automatically stores up to 99 measurement values for each of the 2 users.

Select either user 1 or 2 by pressing the user button .

### Viewing the stored values

Press the M-button  briefly, when the device is switched off. The display first shows «M» , and «A» which stands for the average of all stored values.

Pressing the M-button again displays the previous value. Pressing the M-button repeatedly enables you to move from one stored value to another.

 Blood pressure readings with suboptimal cuff fit -A are not considered in the average value.

 Pay attention that the maximum memory capacity of 99 memories is not exceeded. **When the 99 memory is full, the oldest value is automatically overwritten with the 100 value.** Values should be evaluated by a doctor before the memory capacity is reached – otherwise data will be lost.

### Clearing all values

Make sure the correct user is activated.

If you are sure that you want to permanently remove all stored values, hold down the M-button (the device must have been switched off beforehand) until «CL ALL» appears and then release the button. To permanently clear the memory, press the MAM button while «CL ALL» is flashing. **Individual values cannot be cleared.**

 **Cancel deletion:** press ON/OFF button  while «CL ALL» is flashing.

## 5. Battery indicator and battery change

### Low battery

When the batteries are approximately  $\frac{3}{4}$  empty the battery symbol  will flash as soon as the device is switched on (partly filled battery displayed). Although the device will continue to measure reliably, you should obtain replacement batteries.

### Flat battery - replacement

When the batteries are flat, the battery symbol  will flash as soon as the device is switched on (flat battery displayed). You cannot take any further measurements and must replace the batteries.

1. Open the battery compartment  at the back of the device.
2. Replace the batteries – ensure correct polarity as shown by the symbols in the compartment.
3. To set date and time, follow the procedure described in Section «1. Using the device for the first time».

 The memory retains all values although date and time must be reset – the year number therefore flashes automatically after the batteries are replaced.

### Which batteries and which procedure?

 Use 4 new, long-life 1.5 V, size AA alkaline batteries.

 Do not use batteries beyond their date of expiry.

 Remove batteries if the device is not going to be used for a prolonged period.

## Using rechargeable batteries

You can also operate this device using rechargeable batteries.

-  Only use «NiMH» type reusable batteries.
-  Batteries must be removed and recharged when the flat battery symbol appears. They should not remain inside the device as they may become damaged (total discharge as a result of low use of the device, even when switched off).
-  Always remove the rechargeable batteries if you do not intend to use the device for a week or more.
-  Batteries cannot be charged in the blood pressure monitor. Recharge batteries in an external charger and observe the information regarding charging, care and durability.

## 6. Using a mains adapter

You can operate this device using the Microlife mains adapter (DC 6V, 600 mA).

-  Only use the Microlife mains adapter available as an original accessory appropriate for your supply voltage.
-  Ensure that neither the mains adapter nor the cable are damaged.

1. Plug the adapter cable into the mains adapter socket (7) in the blood pressure monitor.
2. Plug the adapter plug into the wall socket.

When the mains adapter is connected, no battery current is consumed.

## 7. PC-Link functions

This device can be used in conjunction with a personal computer (PC) running the Microlife Blood Pressure Analyzer+ (BPA+) software. The memory data can be transferred to the PC by connecting the monitor via a cable.

If no download-voucher and cable is included download the BPA+ software from [www.microlife.com/software](http://www.microlife.com/software) and use a Micro-USB cable.

## 8. Error Messages

If an error occurs during the measurement, the measurement is interrupted and an error message, e.g. «Err 3», is displayed.

Error	Description	Potential cause and remedy
«Err 1» ①9	Signal too weak	The pulse signals on the cuff are too weak. Re-position the cuff and repeat the measurement.*
«Err 2» ①8-B	Error signal	During the measurement, error signals were detected by the cuff, caused for instance by movement or muscle tension. Repeat the measurement, keeping your arm still.
«Err 3» ①8-C	Abnormal cuff pressure	An adequate pressure cannot be generated in the cuff. A leak may have occurred. Check that the cuff is correctly connected and is not too loose. Replace the batteries if necessary. Repeat the measurement.
«Err 5»	Abnormal result	The measuring signals are inaccurate and no result can therefore be displayed. Read through the checklist for taking a reliable measurement and then repeat the measurement.*
«Err 6»	MAM Mode	There were too many errors during the measurement in MAM mode, making it impossible to obtain a final result. Read through the checklist for taking a reliable measurement and then repeat the measurement.*
«HI»	Pulse or cuff pressure too high	The pressure in the cuff is too high (over 299 mmHg) OR the pulse is too high (over 200 beats per minute). Relax for 5 minutes and repeat the measurement.*
«LO»	Pulse too low	The pulse is too low (less than 40 beats per minute). Repeat the measurement.*

\* Please immediately consult your doctor, if this or any other problem occurs repeatedly.

## 9. Safety, care, accuracy test and disposal

### Safety and protection

- Follow instructions for use. This document provides important product operation and safety information regarding this device. Please read this document thoroughly before using the device and keep for future reference.
- This device may only be used for the purposes described in these instructions. The manufacturer cannot be held liable for damage caused by incorrect application.
- This device comprises sensitive components and must be treated with caution. Observe the storage and operating conditions described in the «Technical Specifications» section.
- The cuffs are sensitive and must be handled with care.
- Only pump up the cuff once fitted.
- Do not use this device if you think it is damaged or notice anything unusual.
- Never open this device.
- Read the additional safety information provided within the individual sections of this instruction manual.
- The measurement results given by this device is not a diagnosis. It is not replacing the need for the consultation of a physician, especially if not matching the patient's symptoms. Do not rely on the measurement result only, always consider other potentially occurring symptoms and the patient's feedback. Calling a doctor or an ambulance is advised if needed.

Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed. Be aware of the risk of strangulation in case this device is supplied with cables or tubes.

### Contra-indications

Do not use this device if the patient's condition meets the following contra-indications, to avoid inaccurate measurements or injuries.

- The device is not intended for measuring blood pressure in pediatric patients of age younger than 12 years old (children, infant, or neonates).
- Presence of significant cardiac arrhythmia during measurement may interfere with blood pressure measurement and affect the reliability of blood pressure readings. Consult with your doctor about whether the device is suitable for use in this case.

- The device measures blood pressure using a pressured cuff. If the measuring limb suffers from injuries (for example open wounds) or under conditions or treatments (for example intravenous drip) making it unsuitable for surface contact or pressurization, do not use the device, to avoid worsening of the injuries or conditions.
- Patient motions during measurement may interfere with the measurement process and influence results.
- Avoid taking measurements of patients with conditions, diseases, and susceptible to environment conditions that lead to uncontrollable motions (e.g. trembling or shivering) and inability to communicate clearly (for example children and unconscious patients).
- The device uses oscillometric method to determine blood pressure. The arm being measure should have normal perfusion. The device is not intended to be used on a limb with restricted or impaired blood circulation. If you suffer with perfusion or blood disorders, consult your doctor before using the device.
- Avoid taking measurement on the arm on the side of a mastectomy or lymph node clearance.
- Do not use this device in a moving vehicle (for example in a car or on an aircraft).

### WARNING

Indicates a potentially hazardous situation, which if not avoided, could result in death or serious injury.

- This device may only be used for the intended uses described in this Instructions for Use. The manufacturer cannot be held liable for damage caused by incorrect application.
- Do not change the patient medication and treatment based the result of one or multiple measurements. Treatment and medication changes should be prescribed only by a medical professional.
- Inspect the device, cuff, and other parts for damage. DO NOT USE the device, cuff or parts if they appear damaged or operating abnormally.
- Blood flow of the arm is temporarily interrupted during measurement. Extended interruption of blood flow reduces peripheral circulation and may cause tissue injury. Beware of signs (for example tissue discoloration) of impeded peripheral circulation if taking measurements continuously or for an extended period of time.

- Prolonged exposure to cuff pressure will reduce peripheral perfusion and may lead to injury. Avoid situations of extended cuff pressurization beyond normal measurements. In the case of abnormally long pressurization, abort the measurement or loose the cuff to depressurize the cuff.
- Do not use this device in oxygen rich environment or near flammable gas.
- The device is not water resistant or water proof. Do not spill or immerse the device in water or other liquids.
- Do not disassemble or attempt to service the device, accessory and parts, during use or in storage. Access to the device internal hardware and software is prohibited. Unauthorized access and servicing of the device, during use or in storage, may compromise the safety and performance of the device.
- Keep the device away from children and people incapable of operating the device. Beware of the risks of accidental ingestion of small parts and of strangulation with the cables and tubes of this device and accessories.
- Protect the device and accessories from the following to avoid damaging the device:
  - water, other liquids, and moisture
  - extreme temperatures
  - impacts and vibrations
  - direct sunlight
  - contamination and dust
- This device has 2 user settings. Please ensure it is cleaned and disinfected between users to prevent any cross contamination.
- Stop using this device and cuff and consult with your doctor if you experience skin irritation or discomfort.

### CAUTION

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient, or cause damage to the device or other property.

- The device is intended only for measuring blood pressure at upper arm. Do not measure other sites because the reading does not reflect your blood pressure accurately.
- After a measurement is completed, loosen the cuff and rest for > 5 minutes to restore limb perfusion, before taking another measurement.
- Do not use this device with other medical electrical (ME) equipment simultaneously. This may cause device malfunction or measurement inaccuracies.
- Do not use this device in proximity of high frequency (HF) surgical equipment, magnetic resonance imaging (MRI) equipment, and computerized tomography (CT) scanners. This may cause device malfunction and measurement inaccuracies.
- Use and store the device, cuff and parts in temperature and humidity conditions specified in the «Technical Specifications». Usage and storage of the device, cuff and parts in conditions outside ranges given in the «Technical Specifications» may results in device malfunction and the safety of usage.

### Electromagnetic Compatibility Information

This device is compliant with EN60601-1-2: 2015 Electromagnetic Disturbances standard.

This device is not certified to be used in vicinity of High Frequency (HF) medical equipment.

Do not use this device close to strong electromagnetic fields and portable radio frequency communication devices (for example microwave oven and mobile devices). Keep a minimum distance of 0.3 m from such devices when using this device.

### Device care

Clean the device only with a soft, dry cloth.

### Cleaning the cuff

The cuff delivered with this device is washable.

1. Remove the cuff connector (11) from the cuff tube (12) and carefully pull the bladder through the opening at the edge of the cuff cover.
2. Hand wash the cuff cover in soapsuds: not hotter than 30 °C.
3. Completely dry the cuff cover by linen drying.
4. Loop the cuff tube back through its opening and carefully place the bladder flat in the cuff cover.
5. Reattach the cuff connector on the cuff tube.

The bladder must lay straight in the cuff cover, not folded.

Do not use fabric softener.

**WARNING:** Do not wash the cuff in a washing machine or dishwasher!

**WARNING:** Do not dry the cuff cover in a tumble dryer!

**WARNING:** Under no circumstances may you wash the inner bladder!

## Accuracy test

We recommend this device is tested for accuracy every 2 years or after mechanical impact (e.g. being dropped). Please contact your local Microlife-Service to arrange the test (see foreword).

## Disposal

Batteries and electronic devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

## 10. Guarantee

This device is covered by a **5 year guarantee** from the date of purchase. During this guarantee period, at our discretion, Microlife will repair or replace the defective product free of charge.

Opening or altering the device invalidates the guarantee.

The following items are excluded from the guarantee:

- Transport costs and risks of transport.
- Damage caused by incorrect application or non-compliance with the instructions for use.
- Damage caused by leaking batteries.
- Damage caused by accident or misuse.
- Packaging/storage material and instructions for use.
- Regular checks and maintenance (calibration).
- Accessories and wearing parts: Batteries, power adapter (optional).

The cuff is covered by a functional guarantee (bladder tightness) for 2 years.

Should guarantee service be required, please contact the dealer from where the product was purchased, or your local Microlife service. You may contact your local Microlife service through our website: [www.microlife.com/support](http://www.microlife.com/support)

Compensation is limited to the value of the product. The guarantee will be granted if the complete product is returned with the original invoice. Repair or replacement within guarantee does not prolong or renew the guarantee period. The legal claims and rights of consumers are not limited by this guarantee.

## 11. Technical Specifications

<b>Operating conditions:</b>	10 - 40 °C / 50 - 104 °F 15-90 % relative maximum humidity
<b>Storage conditions:</b>	-20 - +55 °C / -4 - +131 °F 15-90 % relative maximum humidity

<b>Weight:</b>	402 g (including batteries)
<b>Dimensions:</b>	138 x 94.5 x 62.5 mm
<b>Cuff size:</b>	from 17-52 cm according to the cuff sizes (see «Selecting the correct cuff»)

**Measuring procedure:** oscillometric, corresponding to Korotkoff method: Phase I systolic, Phase V diastolic

**Measurement range:** SYS: 60 - 255 mmHg  
DIA: 40 - 200mmHg  
Pulse: 40 - 199 beats per minute

**Cuff pressure display range:** 0 - 299 mmHg

**Resolution:** 1 mmHg

**Static accuracy:** within  $\pm 3$  mmHg

**Pulse accuracy:**  $\pm 5$  % of the readout value

**Voltage source:**

- 4 x 1.5 V alkaline batteries; size AA
- Mains adapter DC 6V, 600 mA (optional)

**Battery lifetime:** approx. 580 measurements (using new batteries)

**IP Class:** IP 20

**Reference to standards:** IEC 80601-2-30; IEC 60601-1; IEC 60601-1-2 (EMC); IEC 60601-1-11

**Expected service life:** Device: 5 years or 10000 measurements, whichever comes first  
Accessories: 2 years or 5000 measurements, whichever comes first

This device complies with the requirements of the Medical Device Directive 93/42/EEC.

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